


P-5-22 Bed/Floor Sensor Pad 24" x 30" Mounts on mattress or floor.

USER MANUAL QUICK START GUIDE

Components:

1. Sensor Pad with 4 velcro tabs for mounting.

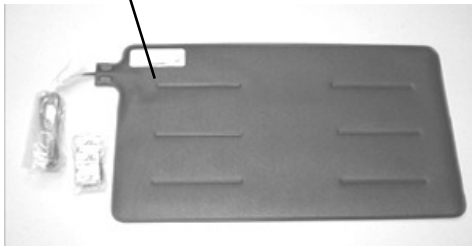

Step 1.

Uncoil the cord from the sensor pad. Place the sensor pad on the bed with the cord coming straight out and down under the mattress, protruding towards the back of the bed towards the alarm unit. The pad would normally be placed on the vinyl mattress just below waist level and the alarm would sound when the patient rose from the bed (weight off). For a more immediate response place the pad under the patients shoulder blades. Then with the patient first attempt to rise (releasing pressure on sensor pad) the alarm will sound.

The RJ22 connection at the end of the cord may be connected into any of Emergency Caller Products alarm units such as the Patient Alert (PA), the Extended Range Transmitter (TR-300), or the Long Range Transmitter (LT-900). Consult your Emergency Caller Products sales representative to learn how these other products work with this model sensor pad. The sensor pad may also be connected into other products that will receive the RJ22 connector.

Sensor Pads detect both Applied Weight "Patient ON" pad or Removed Weight "Patient Off". Selection made by switch on Universal Transmitter Unit or other monitoring products. Can be used to detect patient in bathtub or shower for entry or exit (weight off or weight on). Sensor pads are slip resistant and completely waterproof. Complete with velcro tabs and 8 ft. leads (except Chair-seat Pad with a 4 ft.lead). Lifetime sensor pads replace all disposable sensor pads. Emergency Caller Sensor pads have been used for over 25 years in security applications. They are durable, safe and withstand both outdoor and indoor environments. With life-time guarantee you never again need to be concerned about constant replacement and expensive costs.